Supplement of Aerosol Research, 3, 557–568, 2025 https://doi.org/10.5194/ar-3-557-2025-supplement © Author(s) 2025. CC BY 4.0 License.





# Supplement of

# A simple, versatile approach for coupling a liquid chromatograph and chemical ionization mass spectrometer for offline analysis of organic aerosol

Andre F. Schaum et al.

Correspondence to: Paul J. Ziemann (pazi4520@colorado.edu)

The copyright of individual parts of the supplement might differ from the article licence.

### S1. Sampling and Filter Extraction Methods

Chamber-Generated SOA Filter Extraction

Table S1. Masses of PTFE filters pre- and post-sampling from chamber. All mass acquired was assumed to be SOA formed through the reaction of  $\alpha$ -pinene and O<sub>3</sub>.

Filter No.	Pre-Sampling Mass (mg)	Post-Sampling Mass (mg)	SOA Mass (mg)
1	89.388	90.075	0.687
2	90.021	90.715	0.694

#### OA Sample Collection

Ambient OA was collected at the Denver La Casa CDPHE Air Quality Monitoring Site (4545 Navajo St. Denver, CO) from October 10-18, 2024. The sampler was operated from 04:00-10:00 each day for a total of 42 hours over the study period. Aerosol samples were collected on Pall Tissuquartz filters (P.N. 7204) using an MCV high volume sampler affixed with a PM1 sampling head. The sampler flow rate was set to 30 m<sup>3</sup> h<sup>-1</sup>. At the end of the sampling period, the filters were transported on ice and stored at -10 °C until extraction.

#### Field-Collected OA Filter Extraction

A circular subsection ( $\sim$ 31 cm²) was cut from the filter using a clean razor blade. The filter was then extracted using 12 mL of ethyl acetate over 2.5 days and the resulting solutions were filtered through Teflon filters (1 µm pore size). The filtered extract was dried under a gentle stream of dry N<sub>2</sub> in a preweighed HPLC vial and the final mass was recorded. The sample was reconstituted in HPLC-grade acetonitrile (Fisher Chemical), to the desired concentration and stored at -10 °C until analysis by HPLC-VIA-CIMS.

## S2. Separation and Detection of Carboxylic Acids

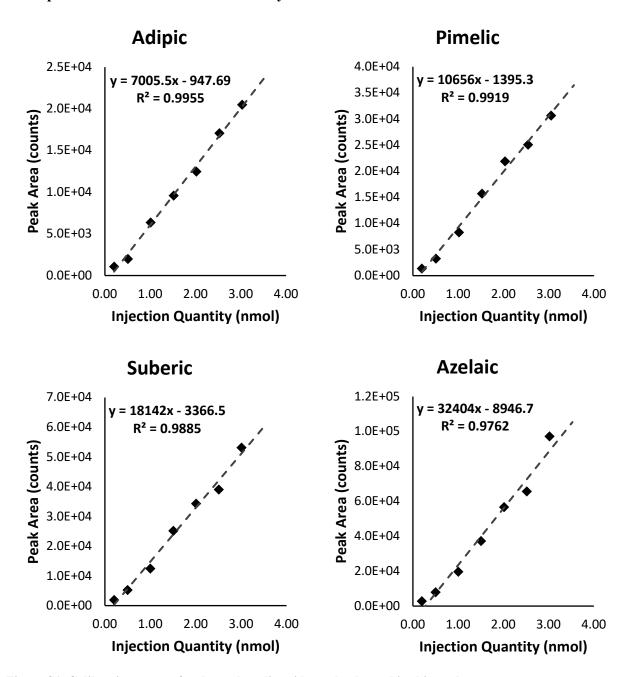


Figure S1. Calibration curves for the carboxylic acid standards used in this study.

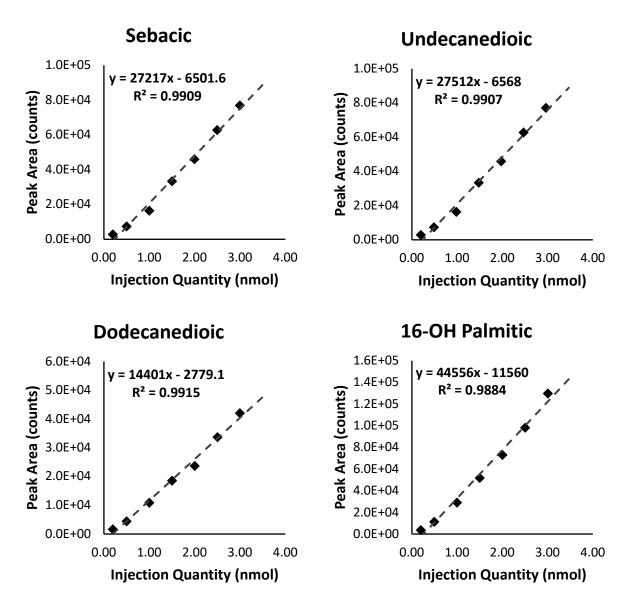


Figure S1 cont. Calibration curves for the carboxylic acid standards used in this study.

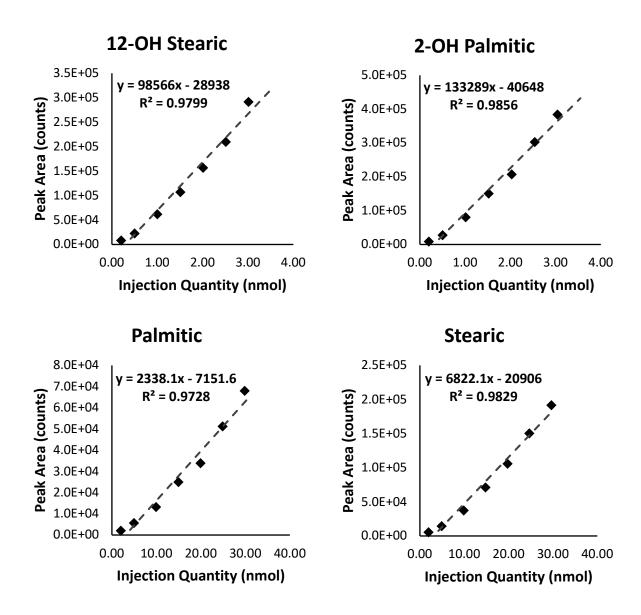


Figure S1 cont. Calibration curves for the carboxylic acid standards used in this study.

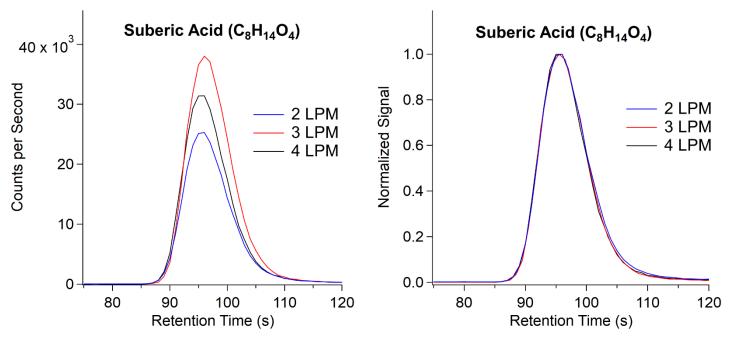


Figure S2. Comparison of raw (left) and normalized (right) suberic acid signals at different system flow rates.

# S3. α-Pinene Ozonolysis SOA

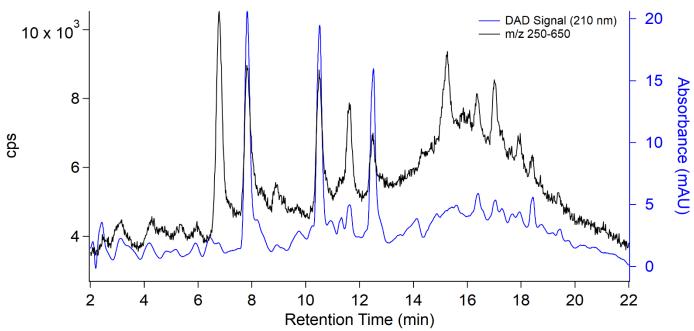


Figure S3. HPLC-VIA-CIMS and DAD chromatograms of  $\alpha$ -pinene ozonolysis SOA. The sum of CIMS signal from m/z 250-650 was used to exclude reagent ion signal while including all expected analyte signals.

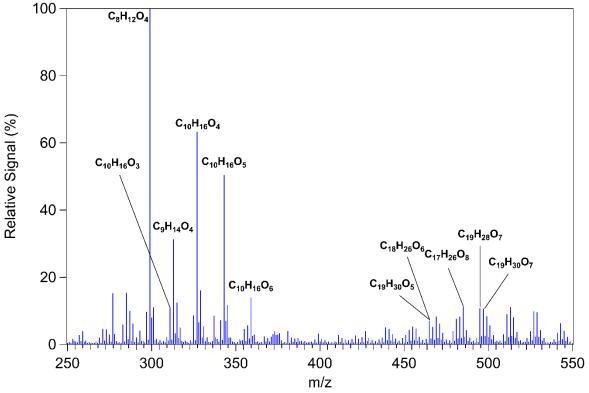


Figure S4. Average mass spectrum of  $\alpha$ -pinene ozonolysis SOA between retention times of  $\sim$ 6 and  $\sim$ 18 min during HPLC-VIA-CIMS analysis. Compounds discussed in the main text are annotated without I<sup>-</sup>. Signals were normalized relative to m/z 299.

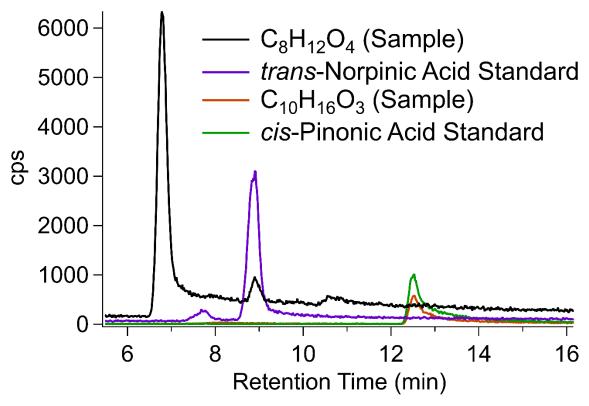


Figure S5. Chromatograms of trans-norpinic ( $C_8H_{12}O_4$ ) and cis-pinonic acid ( $C_{10}H_{16}O_3$ ) standards compared with corresponding chromatograms of  $\alpha$ -pinene SOA. The first peak in the trans-norpinic acid trace (~7.5 min) is likely an impurity with the same molecular formula. The  $C_8H_{12}O_4$  peak at 7 min in the SOA sample is assumed to be a structural isomer of trans-norpinic acid.

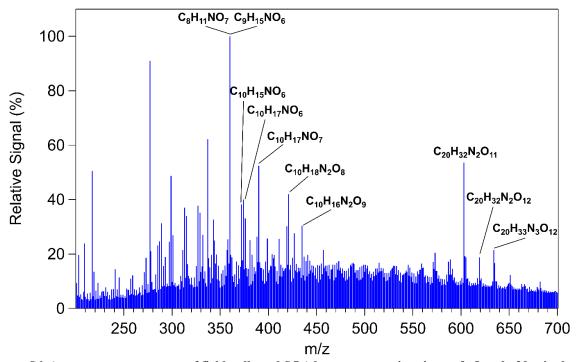


Figure S6. Average mass spectrum of field-collected SOA between retention times of  $\sim$ 8 and  $\sim$ 20 min during HPLC-VIA-CIMS analysis. Compounds discussed in the main text are annotated without I<sup>-</sup>. Signals were normalized relative to m/z 360.

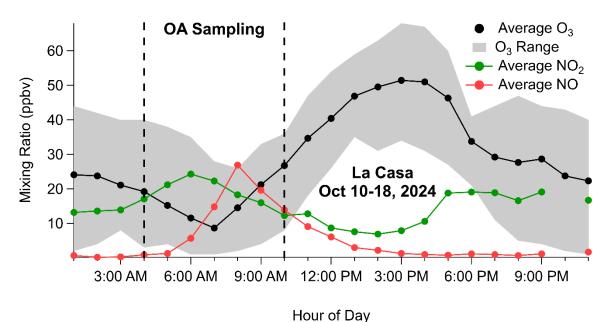


Figure S7. Average hourly concentrations of O<sub>3</sub>, NO<sub>2</sub>, and NO at the Denver, CO ASCENT site during the week of October 10-18, 2024. All data were obtained from the Colorado Department of Public Health & Environment through publicly available air quality reports (https://www.colorado.gov/airquality/report. aspx).