

A rapid semi-quantitative screening method to assess chemicals present in heated e-liquids and e-cigarette aerosols

Natalie Anderson^{1,2}, Paul Pringle³, Ryan Mead-Hunter¹, Benjamin Mullins¹, Alexander Larcombe^{1,2}, Sebastien Allard^{3,4}

5

¹School of Population Health, Curtin University, Perth, Australia, ²Wal-yan Respiratory Research Centre, Telethon Kids Institute, Nedlands, 6009, Australia, ³Curtin Water Quality Research Centre, School of Molecular and Life Sciences, Curtin University, Perth, Australia, ⁴Commonwealth Scientific and Industrial Research Organisation (CSIRO), Land and Water, Floreat, 6014, Perth, Western Australia, Australia.

10

Correspondence to: Natalie Anderson (Natalie.anderson@telethonkids.org.au)

Abstract. Introduction: Electronic cigarettes (e-cigarettes) lack regulatory status as therapeutic products in all jurisdictions worldwide. They are potentially unsafe consumer products, with significant evidence they pose a risk to human health. Therefore, developing rapid, economical test methods to assess the chemical composition of e-liquids in heated and unheated forms and the aerosols produced by e-cigarettes is crucial. **Methods:** Four different e-liquids were heated using two different methods: 1) “typical” vaping using an e-cigarette device, by cycling “on” for three seconds every minute for two hours (e-liquid obtained from remainder in the tank and aerosol collected in an impinger) and, 2) “accelerated” heating, using an e-cigarette coil, submerged in e-liquid, and heating in short 20 second bursts “on” then 20 seconds “off” for two minutes only (liquid traps aerosol produced). All e-liquids were then analyzed to test for the presence and quantity of 13 chemicals by gas-chromatography mass-spectrometry and compared to an unheated sample. **Results:** E-liquids heated with the “accelerated” method showed a comparable trend to the “typical” heating method, i.e. increase or decrease in chemical compound quantity, for more than two-thirds of the detected compounds analyzed over all e-liquids. Six chemicals were detected as aerosol from the impinger fluid with the “typical” heating method, at negligible levels. **Conclusion:** We propose that this accelerated version of the typical vaping method could form the basis of a standardized screening tool to test heated e-liquids (and e-cigarette aerosols) for harmful or banned substances. This will ensure that only approved products reach the consumer and reduce potential e-cigarette harms.

25

1 Introduction

Chemicals are present in the aerosol produced by electronic cigarettes (e-cigarette) when the e-liquid it contains is heated and aerosolised (Goniewicz et al., 2014, European Parliament, 2014). The e-cigarette is referring to the aerosol generating device, which uses the “e-liquid” to create aerosol by evaporation condensation method. If in sufficient quantities, and dependent on the hazard type, chemicals present in an e-cigarette-produced aerosol have the potential to negatively impact health when inhaled (European Association for the Co-ordination of Consumer Representation in Standardisation, 2021). For example,

30

group one carcinogens acetaldehyde and formaldehyde have been found in e-cigarette aerosols and are degradation products of the e-liquid base components, propylene glycol and glycerol (Goniewicz et al., 2014). Other chemical ingredients found in unheated e-liquids that pose a risk to human health (e.g. respiratory irritants, sensitizers etc.) may include flavors or solvents added to the e-liquids. These chemicals can be present in unheated e-liquid, heated e-liquid, aerosolised e-liquid, or any combination of the three forms. Yet, as recently as 2015, no country in the world regulated e-liquid ingredients beyond nicotine levels.

As more evidence emerges suggesting that e-cigarette aerosols negatively impact health, regulation is rapidly evolving in this area (European Parliament, 2014, Budzyńska, Sielemann, Puton, & Surminski, 2020, Therapeutic Goods Administration, 2021). For example, the European Union (EU) Tobacco Product Directive (TPD) states that only ingredients in nicotine-containing e-liquids that do not pose a risk to human health in heated or unheated forms can be used (European Parliament, 2014) and many countries (including United Kingdom, Germany, and France) banned some ingredients as a result (European Parliament, 2014, European Association for the Co-ordination of Consumer Representation in Standardisation, 2021, Budzyńska, Sielemann, Puton, & Surminski, 2020). A modest ingredient ban for nicotine-containing e-liquids came into effect in October 2021 in Australia, prior to a complete ban on non-prescription e-cigarette (and e-liquid) importation and sale in Australia being announced in May 2023 (Nogrady, 2023, Therapeutic Goods Administration, 2021). Known health effects of banned ingredients include e-cigarette or Vaping Associated Pulmonary (or Lung) Injury (aka VAPI, or VALI/EVALI). There is emerging evidence of self-reported lung conditions associated with e-cigarette use (Greenhalgh, 2019, Osei et al., 2020, Bircan, Bezirhan, Porter, Fagan, & Orloff, 2021). Health concerns are amplified by a multibillion-dollar market driving product sales and increasing use in young people. Consequently, pre-market approvals (including unheated ingredient listing) are preferred by regulators which allow independent determination of product safety (Australian Government, 2023, United States Food and Drug Administration, 2023).

Whilst chemical testing of unheated e-liquids is relatively common, the safety of e-liquids remains largely unassessed due to the sheer scale of the market. Additionally, except for the EU TPD guidance (to the best of our knowledge), testing heated e-liquids for inclusion in pre-market product approval is not required in Australia or elsewhere (European Parliament, 2014, Greenhalgh, 2019, Scientific Committee on Health Environmental and Emerging Risks 2021, Larcombe et al., 2021). Methods to assess chemical content generally involve testing of unheated e-liquids or the e-cigarette aerosol generated, but not the heated e-liquid (Scientific Committee on Health Environmental and Emerging Risks, 2021). However, the heated e-liquid is more representative of what the user inhales compared to the unheated e-liquid, and easier to assess when compared to e-cigarette aerosol (Larcombe et al., 2021, Erythropel et al., 2019). Importantly, chemicals present in e-liquids are known to degrade due to heating, either by boiling or evaporative-convection (depending on wetted-wick temperature) and this can be exacerbated by presence of catalytic surfaces such as Kanthal. The secondary products formed, at high (>200°C) or low temperatures (<200 °C), may have increased or decreased toxicity compared to the parent compound (Erythropel et al., 2019, Floyd, Queimado, Wang, Regens, & Johnson, 2019, Goniewicz et al., 2014, Jaegers, Hu, Weber, & Hu, 2021, Li et al., 2021, Zhao, Shu, Guo, & Zhu, 2016, Saliba et al., 2018).

Many approaches have been trialled to address the difficulties in directly assessing the chemical composition of e-cigarette produced aerosol. Such challenges include collecting enough aerosol to perform an assay (and therefore detect potentially toxic compounds), and overload of the main excipients (propylene glycol and vegetable glycerin). These difficulties induce increasingly complicated test methods, yet current methods have been shown to be outdated or limited (Floyd, Queimado, Wang, Regens, & Johnson, 2019, Herrington & Myers, 2015, Scientific Committee on Health Environmental and Emerging Risks, 2021). Despite the expansive range of tests, a simple test capable of assessing heated and unheated e-liquids, and the aerosol produced all at once, is yet to be established, but vital.

In this study we aimed to: 1) establish and validate a simple “accelerated”, method of heating and aerosolising e-liquids (i.e. “vaping”) that would be comparable to, but quicker than, the “typical” method of heating and aerosolising e-liquid and; 2) assess the chemical composition of the heated e-liquid and the aerosol produced all at once. We hypothesized that the “accelerated” method and the commonly used “typical” vaping method would result in similar heating-induced chemical changes in the e-liquids. The two heating method were compared by measuring the presence and concentration of a range of chemicals. We hope that this “accelerated” test methodology can form the critical first step in establishing a rapid test for screening of e-liquids for banned substances.

2 Materials and Methods

2.1 “Typical” vaping process

A set up was designed to replicate the heating/cooling process an e-liquid would undergo when an e-cigarette is used in a “typical” way (Figure 1 (a)) (Etter & Bullen, 2014, St Helen et al., 2016, Cooperation Center for Scientific Research Relative to Tobacco (CORESTA), (2015)). The method allowed sample collection at two points for analysis of heated e-liquid and aerosol respectively: 1) from the remainder in the e-liquid tank (atomiser), and 2) from the impinger (Figure 1 (a)). To begin, the e-cigarette (MVP4, Innokin, Shenzhen, China, operating wattage range 6–100 W, temperature range 150–315 °C, maximum current 35.5A) atomizer (the e-liquid tank, containing e-cigarette coil) was filled with ~3.5 mL of e-liquid and the impinger was filled with 5 mL of e-liquid excipient (50:50 glycerin–propylene glycol (v/v)), Sigma Aldrich, Milwaukee WI, USA). A flow of ~3 L/min ambient filtered air was drawn through the system via laboratory bench vacuum and kept stable through monitoring with a flow meter (Max 5 Lpm, TSI, 800669, Shoreview, MN). New coils (Kanthal BVC, 100–200 W, 0.28 Ω, Innokin Scion) were used each time to avoid cross contamination of chemical species and to control for coil ageing effects. The e-cigarette device was set to 80 W (reading 0.28–0.35 Ω) each time the device was connected to the atomiser. To vape the device in line with recommendations by the Cooperation Centre for Scientific Research Relative to Tobacco (CORESTA), the ignition button was held for ~3 seconds, and the aerosol drawn from the device into a 60 mL syringe and then expelled through two ~ 4 mm ID, 15 cm tubing lengths into the custom made (27 L – 30x30x30 cm) chamber, using a three-way tap (Dispoflex™, Disposafe health and life care Ltd, Haryana, India) (Cooperation Center for Scientific Research Relative to Tobacco (CORESTA) (2015)). This process was repeated every minute for two hours (with the atomizer tank

refilled after ~60 minutes), for 120 puffs total. While we acknowledge that vaping topography is extremely variable, 120 puffs
100 over a 2 hour period (120 x 60 mL, puffs, therefore 7.2 L of inhaled aerosol containing air) was chosen to be representative of
what a typical vaper might use (Etter, 2014). After heating, the liquid was transferred to glass sample vials, and kept at 4°C to
minimise the loss of volatiles.

2.2 “Accelerated” vaping process

An accelerated ageing/vaping process was developed, based on standard tests for ageing/oxidation of oils (Figure 1(b))
105 (American Society for Testing Materials, 2009). Our premise for collection of e-cigarette aerosols in the liquid was as follows:

1. An e-cigarette is an evaporation condensation aerosol generator – intended to modify the e-liquid as little as possible during aerosolization, however, it does thermo-oxidise, hence the need for this research;
2. Our “accelerated” method of heating the e-liquid via a submerged coil creates a “bubbling aerosol generator” (Vidamantas, 1997). Like an evaporation aerosol condensation generator, a bubbling generator will modify the e-
110 liquid minimally, however, may allow more volatile compounds to preferentially aerosolise;
3. The creation of an aerosol via bubbling can allow aerosol capture either whilst bubbling through the bulk liquid (when cooling) or at the gas-liquid surface (Ghiaassiaan, 1997 and Koch, 2012);
4. Surface bubbles can generate aerosol either by jet or film droplets when they burst, and based on combinations of surface tension and bubble size, aerosol will recombine with the liquid the bubble arises from when it bursts (Koch,
115 2012, Mead-Hunter, 2018).
5. Thereby, through a combination of these processes it is reasonable to assume we retain a representative sample of the same material that is aerosolised, as well as possibly more of the thermo-oxidised (aged) material we are interested in.

To create the bubbler, an identical Kanthal BVC coil, as used in the “typical” vaping process, was connected to a power supply
120 (MP3090, PowerTech, China), by means of solid copper wires connected to each end of the coil (end cap removed). The power supply was set at 7.4 V and 27 A (0.274 Ω) to stay within the maximum power (200 W) of the coil used for the “typical” process and to ensure that the resistance matched that of the “typical” vaping process. The coil was placed in 100 mL beaker, which was open to air, held on a 45° angle with a clamp-stand and ~30 mL of e-liquid was poured into the beaker, enough to completely submerge the coil and ensure the full volume of liquid would not heat to boiling temperature within the one minute
125 total heating period (20 seconds on 20 seconds off x 3), and limited planar surface evaporation would occur (Figure 1 (b)). The 45° angle was used both to minimize the liquid volume needed to immerse the coil and to ensure any aerosol (or vapor) produced would recondense on the wall of the beaker allowing it to be collected for sampling. The power supply was then turned on to operate the coil for 5 x 20 second “burst” intervals with 20 second pauses interspersed for a total “on” time of 1 minute, mimicking a “short cluster” vaping pattern for a user (St Helen et al., 2016). After heating, the liquid was transferred
130 to glass sample vials, and kept at 4°C to minimise the loss of volatiles.

2.3 Sample and chemical selection

Four flavoured e-liquids, labelled “nicotine-free”, were assessed – “Butterscotch Tobacco”, “Menthol”, “Choc Caramel”, and “Tiramisu” which were purchased from online suppliers and analyzed as 50:50 propylene glycol–glycerin (v/v) ratios. The

propylene glycol: glycerol mixture was selected as 50:50 since it is a commonly sold ratio. Each e-liquid chemical composition was assessed using both methods to quantify 13 chemicals: 4-(4-methoxyphenyl)-2-butanone, ethyl vanillin, eugenol, nicotine, menthol, thymol, ethyl maltol, trans-cinnamaldehyde, 2-chloro-phenol, benzyl alcohol, benzaldehyde, and furfural, with a molecular weight range from 178.23 to 96.09 g/mol. The 13 chemicals were chosen based on (i) being previously identified, known ingredients in e-liquids (ii) the availability of a standard for the chemical (Larcombe et al., 2021).

2.4 Chemical analysis method of “Accelerated” and “Typical” vaping process.

Thirteen chemicals were tested for, in four different e-liquids, using gas-chromatography mass-spectrometry. For each of the four e-liquids, we tested for chemicals in three forms – (i) “unheated” e-liquid (i.e. straight out of the bottle), (ii) remainder of e-liquid in the atomizer and collected from the impinger after “typical” vaping and (iii) e-liquid remaining in the beaker after “accelerated” vaping. The latter sample (iii) was taken in order to detect aerosols, and assuming that aerosols (not vapor) would be captured in the e-liquid with/during the accelerated method. The aerosol generated from the “typical” method was captured in an impinger containing 50:50 (v/v) glycerin–propylene glycol. Our intention was for this collected aerosol in the impinger to be added to the atomiser tank sample, for equivalent comparison to accelerated sample, however negligible values for the impinger result meant that they were excluded from the final analysis.

Samples obtained from both methods used to heat e-liquids were compared to unheated e-liquids, both within e-liquid type, and within chemical compound, with the purpose of the comparison being to identify trends of increase or decrease from unheated e-liquid.

2.5 Chemical detection and analysis

Chemical analysis of “accelerated” and “typical” vaping process e-liquids, including sample and chemical detection, has been previously described in detail elsewhere (Larcombe et al., 2021). The samples (0.25 g) were accurately weighed and placed into amber vials with 4.75 mL ultrapure water. Thereafter, 10 μ L of a 1 g/L 4-bromophenol-d4 stock solution was added as an internal standard. Prior to the analysis, 1.6 g of analytical grade sodium chloride was added to increase volatilisation and the vials tightly capped. To facilitate adsorption, the samples were incubated at 90°C for 15 min prior to solid-phase micro-extraction using a divinylbenzene/carboxen/polydimethylsiloxane fiber from Supelco® allowing for 13.6 min adsorption of the analytes on the fiber. The fiber was then desorbed at 250 °C in the injector in spitless mode for 5 min followed by 15 min in split mode. A Gerstel MPS2 multifunction autosampler was used to perform automated solid-phase micro-extraction injections. Analysis were carried out with an Agilent 6890N gas-chromatograph interfaced with an Agilent 5973 Network Mass Selective Detector, fitted with a HP-INNOWax polyethylene glycol stationary phase capillary column (30 m; 0.25 mm; 0.25 μ m, Agilent J&W, Australia), to separate polar compounds. A constant flow (1.2 mL.min⁻¹) of helium (99.999% pure, BGC, Australia) was used as a carrier gas. Optimal gas-chromatography mass-spectrometry conditions were determined, as measured by maximum sensitivity, baseline separation of analytes and gaussian peak shapes. In order to ensure a good separation of the different analytes, the oven was held isothermal at 37 °C (2 min), then heated to 260 °C at 5 °C.min⁻¹, and

held at the final temperature for 10 min. Detection of analytes was carried out using a mass spectrometer in electron impact ionization mode at 70 eV. The mass spectrometer quadrupole temperature was set at 150°C and the mass spectrometer source at 230°C. The compounds were identified using a combination of their retention times, comparison of the mass spectra data of pure compounds and the specific diagnostic ion fragments of each component, with the National Institute of Standards and
170 Technology Mass Spectral search program from the NIST/EPA/NIH EI and NIST Tandem Spectral Library which came integrated with the analysis software.

3 Results

Over all e-liquids, three of the thirteen compounds tested for were not detected in any e-liquid form (4-(4-methoxyphenyl)-2-butanone, thymol, 2-chlorophenol) (Table 1).

175 3.1 Inclusion and exclusion criteria for analysis

There were 16 instances where a chemical was detected in unheated and both heated forms (Table 1) and all were included in analysis. Analysis involved: 1) simple comparison in table format of the heated (two methods) and unheated form of an e-liquid sample and, 2) comparison via fold change compared to unheated ($(Y-X)/X$, where X is the unheated sample (mg/L concentration) and Y is the heated sample (mg/L concentration)) for both “typical” and “accelerated” heating methods.

180 Fold-change analyses were not possible on the following: one chemical was undetected in unheated form but detected in both heated forms (menthol in butterscotch tobacco); one chemical was undetected in unheated form and detected in only one heated form (trans-cinnamaldehyde in choc caramel). A further three chemicals were detected in unheated form and only one heated form (benzyl alcohol (Tiramisu), eugenol (Tiramisu) and furfural (Tiramisu)) (Table 1, represented by italicised values).

3.2 Behaviour of the different chemicals detected in e-liquids

185 To compare the effect of heating, results are displayed as fold change compared to unheated ($(Y-X)/X$, where X is the unheated sample (mg/L concentration) and Y is the heated sample (mg/L concentration)) for both “typical” and “accelerated” heating methods (Figure 2). Specific chemicals (benzaldehyde, benzyl alcohol, ethyl vanillin, ethyl maltol, furfural, menthol, nicotine, and nicotine) were present in unheated form, and both heated forms in 16 instances.

Over all e-liquids, in these 16 instances, 70% (11/16) demonstrated a consistent trend within chemical type i.e. both methods
190 of heating either increased or decreased in concentration compared to unheated sample. Ethyl vanillin (choc caramel and tiramisu), furfural (butterscotch tobacco), ethyl maltol (tiramisu), and benzaldehyde (choc caramel) are the five exceptions.

3.3 Chemical characterization by e-liquid type

In the “Menthol” e-liquid, of the 13 chemicals tested, nine were not detected in the heated or unheated sample (Table 1). Of the four that were detected (nicotine, nicotine, menthol and benzaldehyde), all (4/4, 100%) exhibited the same trend

195 (increasing concentration) after heating when compared to the unheated sample for both “typical” and “accelerated” heating methods (Table 1, Figure 2).

In the “Butterscotch Tobacco” e-liquid, of 13 chemicals tested, eight were not detected in heated or unheated form, and an additional one (menthol) was undetected in unheated form (Table 1). Of the four chemicals detected in each sample, three (3/4, 75%) (ethyl vanillin, benzyl alcohol and benzaldehyde) exhibited the same trend for both “typical” and “accelerated” heating
200 methods when compared with the unheated e-liquid (Table 1, Figure 2). However, furfural increased after “typical” heating but decreased with the “accelerated” method when compared to the unheated sample.

In the “Tiramisu” e-liquid, of 13 chemicals tested, seven were not detected in heated or unheated form, and an additional three were undetected in one form of heating (Table 1). Of the three detected in each sample, one (1/3, 33%) (benzaldehyde) exhibited the same trend after both “typical” and “accelerated” heating (Table 1, Figure 2). The remaining two chemicals
205 detected (ethyl vanillin and ethyl maltol) both decreased after “typical” heating but increased after “accelerated” heating when compared to the unheated sample.

In the “Choc Caramel” e-liquid, of 13 chemicals tested, seven were not detected in heated or unheated form and an additional one (trans-cinnamaldehyde) was undetected in both unheated and heated form (Table 1). Of the five chemicals detected in each sample, three (3/5, 60%) (benzyl alcohol, ethyl maltol, and furfural) exhibited the same trend after both “typical” and
210 “accelerated” heating methods (Table 1, Figure 2). Benzaldehyde increased after “typical” heating but decreased after “accelerated” heating compared to the unheated sample. Conversely, ethyl vanillin increased after “accelerated” heating, but decreased after “typical” heating when compared to the unheated sample.

3.4 Impinger results from “typical vaping” heating method

Only six of the 13 chemicals tested were detected at negligible quantities in the impinger fluid; furfural, benzaldehyde, menthol,
215 benzyl alcohol, ethyl maltol and ethyl vanillin and as such they were unable to be included in final analysis.

4 Discussion

The “accelerated” method used here is simple, cost effective, and has the potential to produce heated e-liquid and aerosol in a similar manner to an e-cigarette but on an accelerated timeframe, allowing chemical assessment in a single experiment. Due to the prohibitive costs of commercially available vaping machines, many “in house” simplified methods and set-ups (e.g. e-
220 cigarette puffing machines) have been developed. However, to the best of our knowledge, they all focus on generation of e-cigarette aerosol and not on assessment of the heated e-liquid (Palazzolo, Caudill, Baron, & Cooper, 2021).

Comparison of the accelerated and typical heated samples with their unheated counterpart showed that in ~70% of the e-liquids tested the heating methods demonstrated a similar trend, i.e. increase or decrease of chemical concentration. The four chemicals implicated in the five differing comparisons were mostly aldehydes (ethyl vanillin (2/11), furfural (1/11), and benzaldehyde
225 (1/11)) except for ethyl maltol (1/11) being an alcohol. In three cases; ethyl vanillin (2/11) and ethyl maltol (1/11), an increase

in chemical concentration with the “accelerated” heating method compared to the “typical” method was detected. The observed “increase” with the “accelerated” method can be attributed to a loss of aerosol with the “typical” vaping experimental method. Rainout (recondensation of the aerosol as it cools) of the liquid aerosol was observed in the three-way tap system and the very thin tubing (ID ~ 4 mm) connecting the tap system to the 27 L chamber (Figure 1, a). A modification of the design (e.g. larger tubing) would reduce the rainout losses, thus increasing the yield from the impinger with the “typical” method and improving reproducibility between methods. It is also possible that the flavor aldehydes were present in their propylene glycol acetal form instead of their aldehyde form, as aldehydes are known to form acetals readily (Erythropel et al., 2019). The inability to fully capture the aerosol from the “typical” heating method due to rainout, meant we are unable to confirm the suitability of our “accelerated” heating method as an impinger for aerosol, but only its validity to compare heating methods.

The remaining two discrepancies involved furfural (1/11) and benzaldehyde (1/11), and these compounds were found in increased quantities with the “typical” method compared to the “accelerated” method. Because these two discrepancies contain low molecular weight/low boiling point products we suspect they may have evaporated more readily (compared to the “typical” method) and that our “accelerated” method was simply unable to capture compounds with low boiling point that volatilise easily (Erythropel et al., 2019). For example, furfural has the lowest boiling point and molecular weight of all chemicals detected (162°C, 96.09 g/mol) and was detected as decreased from the unheated sample in both “butterscotch” and “tiramisu” flavor. While our study design angled the beaker at 45° to allow re-condensation of any vapor on the beaker wall, the experiment was carried out in a ventilated fume hood, which may have increased the loss of highly volatile compounds. In future studies, the addition of a lid on the angled beaker as well as monitoring the liquid temperature on the wick (or surrounding the wick), may help to reduce discrepancies to allow full validation of the method for detection of aerosols (Li et al., 2021, Palazzolo, Caudill, Baron, & Cooper, 2021, Bitzer et al., 2018). Further studies should also be designed to consider the solubility of the flavoring compound in the base excipients (i.e. propylene glycol or glycerol), and not the boiling point, as this is suggested to be a major indicator of whether a compound will be detected in/carried over into an aerosol (Erythropel et al., 2019).

Undetected chemicals included 4-(4-methoxyphenyl)-2-butanone, thymol, and 2-chloro-phenol. Considering no “fruity” flavors were assessed, it is less surprising that 4-(4-methoxyphenyl)-2-butanone was undetected, as it is a raspberry ketone methyl ether – a common flavoring in “berry” flavored e-liquids. Although we were looking for thymol because of its use as a precursor for racemic menthol (produced from m-cresol) its absence may be explained because pulegone and other terpenoids are also used as the precursor (Dylong, Hausoul, Palkovits, & Eisenacher, 2022). Additionally, menthol was not detected in either “tiramisu” or “choc-caramel” flavors so perhaps synthetic analogues were present (where menthol was not detected) such as N-ethyl 2-isopropyl-5-methylcyclohexanecarboxamide (trade name WS-3). As the demand for menthol increases, alternative methods to produce L-menthol are on the rise, such as from citronellal (Dylong, Hausoul, Palkovits, & Eisenacher, 2022).

Chlorophenols like 2-chloro-phenol have previously been detected in e-liquids probably because they are notorious environmental contaminants. Particularly, 2-chloro-phenol is a priority contaminant in both the US and EU and has previously been found in e-liquids (Larcombe et al., 2021, Chivers, Janka, Franklin, Mullins, & Larcombe, 2019, Igbinsosa et al., 2013).

260 Two-chloro-phenol is used for many applications, predominantly as a detergent but also as an intermediate in the
manufacturing of agricultural chemicals, pharmaceuticals, biocides, and dyes. Therefore, it is commonly detected in
environmental water samples after being discharged from industrial effluents (Igbinosa et al., 2013, Yahaya, Okoh, Agunbiade,
& Okoh, 2019). It has previously been suggested that 2-chloro-phenol may be a contaminant from the glycerin excipient, for
two reasons; 1) vegetable glycerin is made from plant crops such as canola and 2-chloro-phenol has been found in canola as a
265 pesticide residue and, 2) glycerin (not from plants) is a by-product of bio-diesel production and biodiesel can be made with
canola (Abdel-Gawad, H., & Hegazi, B., 2010, Yahaya, Okoh, Agunbiade, & Okoh, 2019). It is possible that derivatives of 2-
chloro-phenol or other phenolic derivatives known to be priority contaminants were present. However, it was not within the
scope of the study to assess these compounds as they are not commonly reported to be found/tested for in e-liquids. Eugenol
and trans-cinnamaldehyde were the least detected compounds, they were found in only one flavor e-liquid. Cinnamaldehyde
270 might not be detected because it is known to form propylene glycol acetals ($\alpha\beta$, unsaturated aldehyde) and this form is known
to be relatively reactive (Erythropel et al., 2019).

4.1 Limitations and future directions

A limitation of this study is that we tested a pre-determined list of chemicals, based on our knowledge of known e-liquid
ingredients, available standards, and available analytical methods, rather than looking for a complete chemical characterization.
275 This approach allowed us to test a larger range of e-liquids and demonstrate the utility of the “accelerated” aging technique,
as per the over-arching goal of the study. However, an “open ended” approach may be useful for future studies if this method
is to become standardized. An open-ended approach would allow a more complete comprehension of the aging process and
oxidation reactions occurring but would require a broader range of analytical techniques than demonstrated here. Furthermore,
whilst we assessed for some ingredients which are now banned (i.e. benzaldehyde and cinnamaldehyde), this study was
280 designed and conducted prior to the enactment of banned ingredients by the Therapeutic Goods Administration in 2021
(Therapeutic Goods Administration, 2021) and expanding future analysis to include the full range of banned chemical products
would be helpful. The same comment applies for chemicals banned in other jurisdictions.

The method we describe in this study has many advantages over current methods for testing e-liquids. It is a rapid and
inexpensive set-up allowing assessment of the chemical composition of heated e-liquids and with minor modifications, their
285 resultant aerosols. It could be used with any available coil that can be modified and powered as described. Furthermore, the
accelerated method is likely to capture aerosol generated from a heated e-liquid in a manner comparable to the “typical user”
vaping method as described in CORESTA in terms of both type and quantity of chemicals produced. Our submerged, rapid
heating and cooling method can economically sample heated liquid and aerosols (but not vapor) within a single sample in two
minutes, which may have advantages over some other methods. This method is more representative of what the user inhales
290 as it is testing a heated liquid during exposure to the coil (potentially) catalytic surface, rather than only an unheated e-liquid.

5 Conclusion

In summary, the accelerated method described here is a suitable screening tool for rapid chemical assessment of heated e-liquids and their aerosols, that mimics typical e-cigarette vaping on an accelerated timeframe. The EU TPD recommendation is to assay heated e-liquids, however (to the best of our knowledge) there has only been one previously published study on the effects of aging/heating on e-liquids. We propose that this method (with our recommended improvements) can be used as a standardized screening tool for e-liquids, and their aerosols, to identify potentially harmful chemicals, such as those recently banned in Australia or previously banned in Europe and the United Kingdom. With minor modification, this test could be used prior to importation or sale, to ensure that only tested products, containing approved ingredients, reach the consumer. In the absence of an approved therapeutic goods status for e-cigarettes, the type of high-throughput testing described here is necessary as a minimal precaution to assess and reduce the potential harms of a consumer product that is generally accepted in the public to be a less harmful alternative to smoking.

CRediT statement: Natalie Anderson: data curation, formal analysis, investigation, methodology, project administration, visualization, writing – review and editing, writing – original draft.; Paul Pringle: data curation, formal analysis, investigation, methodology, validation, writing – review and editing.; Ryan Mead-Hunter: formal analysis, investigation, writing – review and editing.; Benjamin Mullins: conceptualization, methodology, resources, writing – review and editing, funding acquisition.; Alexander Larcombe: conceptualization, methodology, project administration, resources, supervision, visualization, writing – review and editing, funding acquisition.; Sebastien Allard: data curation, formal analysis, investigation, methodology, resources, supervision, writing – review and editing, validation, visualization

Acknowledgements: This work was supported by funding from the Minderoo Foundation, Lung Foundation Australia, Cancer Council WA, and the Scottish Masonic Charitable Foundation

Declaration of Interests: All authors declare no competing interest

References

- Abdel-Gawad, H., & Hegazi, B. Fate of ¹⁴C-ethyl prothiofos insecticide in canola seeds and oils. *J ENVIRON SCI HEAL B*, 330 45(2), 116–122. <https://doi.org/10.1080/03601230903471878>, 2010.
- American Society for Testing Materials. Standard test method for bench oxidation of engine oils by ROBO apparatus D7528-09. 2009
- Australian Government, Department of Health and Aged Care, Therapeutic Goods Administration. Webpage content accessed 21 July 2023. Retrieved from: [https://consultations.tga.gov.au/medicines-regulation-division/proposed-reforms-to-the-](https://consultations.tga.gov.au/medicines-regulation-division/proposed-reforms-to-the-regulation-of-nicotine-vap/)
335 [regulation-of-nicotine-vap/](https://consultations.tga.gov.au/medicines-regulation-division/proposed-reforms-to-the-regulation-of-nicotine-vap/)
- Bircan, E., Bezirhan, U., Porter, A., Fagan, P., & Orloff, M. S. Electronic cigarette use and its association with asthma, chronic obstructive pulmonary disease (COPD) and asthma- COPD overlap syndrome among never cigarette smokers. *TOB INDUC DIS*, 19(April), 1-9. doi:10.18332/tid/132833, 2021.
- Budzyńska, E., Sielemann, S., Puton, J., & Surminski, A. L. R. M. Analysis of e-liquids for electronic cigarettes using GC-
340 IMS/MS with headspace sampling. *TALANTA*, 209, 120594. doi:<https://doi.org/10.1016/j.talanta.2019.120594>, 2020.
- Chivers, E., Janka, M., Franklin, P., Mullins, B., & Larcombe, A. Nicotine and other potentially harmful compounds in "nicotine-free" e-cigarette liquids in Australia. *MED J AUSTRALIA*, 210(3), 127-128. doi:10.5694/mja2.12059, 2019.
- Dylong, D., Hausoul, P. J. C., Palkovits, R., & Eisenacher, M. Synthesis of (–)-menthol: Industrial synthesis routes and recent development. *FLAVOUR FRAG J*, 37(4), 195-209. doi:<https://doi.org/10.1002/ffj.3699>, 2022.
- 345 Erythropel, H. C., Jabba, S. V., DeWinter, T. M., Mendizabal, M., Anastas, P. T., Jordt, S. E., & Zimmerman, J. B. Formation of flavorant-propylene Glycol Adducts With Novel Toxicological Properties in Chemically Unstable E-Cigarette Liquids. *NICOTINE TOB RES: official journal of the Society for Research on Nicotine and Tobacco*, 21(9), 1248-1258. doi:10.1093/ntr/nty192, 2019.
- Etter, J.-F., & Bullen, C. A longitudinal study of electronic cigarette users. *ADDICT BEHAV*, 39(2), 491-494, 2014.
- 350 European Association for the Co-ordination of Consumer Representation in Standardisation (ANEC). Position Paper: e-cigarettes and e-liquids - limits for chemicals - basis for discussion. 2021
- European Parliament. Directive 2014/40/EU of the European Parliament and of the Council of 3 April 2014 on the Approximation of the Laws, Regulations and Administrative Provisions of the Member States Concerning the Manufacture, Presentation and Sale of Tobacco and Related Products and Repealing Directive 2001/37/EC (30.12.2018).
- 355 Floyd, E. L., Queimado, L., Wang, J., Regens, J. L., & Johnson, D. L. Electronic cigarette power affects count concentration and particle size distribution of vaping aerosol. *PLOS ONE*, 13(12), e0210147. doi:10.1371/journal.pone.0210147, 2019.
- Ghiaasiaan, S. M., & Yao, G. F. A Theoretical Model for Deposition of Aerosols in Rising Spherical Bubbles due to Diffusion, Convection, and Inertia, *AEROSOL SCI TECH*, 26:2, 141-153, DOI: 10.1080/02786829708965420, 1997.

- Goniewicz, M. L., Knysak, J., Gawron, M., Kosmider, L., Sobczak, A., Kurek, J., . . . Benowitz, N. Levels of selected
360 carcinogens and toxicants in vapour from electronic cigarettes. *TOB CONTROL*, 23(2), 133-139. doi:10.1136/tobaccocontrol-2012-050859, 2014.
- Greenhalgh, E., & Scollo, MM. InDepth 18B: Electronic cigarettes (e-cigarettes). Tobacco in Australia: Facts and issues. Melbourne: Cancer Council Victoria; 2019. Retrieved from <http://www.tobaccoin australia.org.au/chapter-18-harm-reduction/indepth-18b-e-cigarettes>, 2019.
- 365 Herrington, J.S & Myers, C. Electronic cigarette solutions and resultant aerosol profiles, *J CHROMATOGR A*, Volume 1418, Pages 192-199, ISSN 0021-9673, <https://doi.org/10.1016/j.chroma.2015.09.034>, 2015.
- Igbinsola, E. O., Odjadjare, E. E., Chigor, V. N., Igbinsola, I. H., Emoghene, A. O., Ekhaise, F. O., . . . Idemudia, O. G. (2013). Toxicological profile of chlorophenols and their derivatives in the environment: the public health perspective. *SCI WORLD J*, 2013, 460215-460215. doi:10.1155/2013/460215
- 370 Jaegers, N. R., Hu, W., Weber, T. J., & Hu, J. Z. Low-temperature (<200 °C) degradation of electronic nicotine delivery system liquids generates toxic aldehydes. *SCI REP-UK*, 11(1), 7800. doi:10.1038/s41598-021-87044-x, 2021.
- Koch, D & Weber, A. P., Separation of gas-borne nanoparticles in bubble columns. *J AEROSOL SCI.* 53:61-75. doi:10.1016/j.jaerosci.2012.05.012, 2012.
- Larcombe, A., Allard, S., Pringle, P., Mead-Hunter, R., Anderson, N., & Mullins, B. Chemical analysis of fresh and aged
375 Australian e-cigarette liquids. *MED J AUSTRALIA*. doi:10.5694/mja2.51280, 2021.
- Li, Y., Burns, A. E., Tran, L. N., Abellar, K. A., Poindexter, M., Li, X., . . . Nguyen, T. B. Impact of e-Liquid Composition, Coil Temperature, and Puff Topography on the Aerosol Chemistry of Electronic Cigarettes. *CHEM RES TOXICOL*, 34(6), 1640-1654. doi:10.1021/acs.chemrestox.1c00070, 2021.
- Mead-Hunter, R., Gumulya, M.M., King, A. J. C., Mullins, B. J., Ejection of Droplets from a Bursting Bubble on a Free Liquid
380 Surface-A Dimensionless Criterion for "Jet" Droplets. *LANGMUIR*. May;34(21):6307-6313. DOI: 10.1021/acs.langmuir.8b00664. PMID: 29737855, 2018.
- Nogrady B. Australia bans all vapes except on prescription to stem use in children *BMJ OPEN*; 381 :p1014 doi:10.1136/bmj.p1014, 2023.
- Osei, A. D., Mirbolouk, M., Orimoloye, O. A., Dzaye, O., Uddin, S. M. I., Benjamin, E. J., . . . Blaha, M. J. Association
385 Between E-Cigarette Use and Chronic Obstructive Pulmonary Disease by Smoking Status: Behavioral Risk Factor Surveillance System 2016 and 2017. *AM J PREV MED*, 58(3), 336-342. doi:10.1016/j.amepre.2019.10.014, 2020.
- Palazzolo, D. L., Caudill, J., Baron, J., & Cooper, K. Fabrication and Validation of an Economical, Programmable, Dual-Channel, Electronic Cigarette Aerosol Generator. *INT J ENV RES PUB HE*, 18(24), 13190. doi:10.3390/ijerph182413190, 2021.
- 390 Saliba, N. A., El Hellani, A., Honein, E., Salman, R., Talih, S., Zeaiter, J., & Shihadeh, A. Surface Chemistry of Electronic Cigarette Electrical Heating Coils: Effects of Metal Type on Propylene Glycol Thermal Decomposition. *J ANAL APPL PYROL*, 134, 520-525. doi:10.1016/j.jaap.2018.07.019, 2018.

Scientific Committee on Health Environmental and Emerging Risks (SCHEER). Scientific Opinion on electronic cigarettes, 2021.

395 St Helen, G., Ross, K. C., Dempsey, D. A., Havel, C. M., Jacob, P., 3rd, & Benowitz, N. L. Nicotine Delivery and Vaping Behavior During ad Libitum E-cigarette Access. *TOB REGUL SCI*, 2(4), 363-376. doi:10.18001/TRS.2.4.8, 2016.

Therapeutic Goods Administration. (2021). Therapeutic Goods Standard for Nicotine Vaping Products TGO 110 Order. Retrieved from <https://www.legislation.gov.au/Details/F2021L00595>, 2021

400 United States Food and Drug Administration. Webpage content accessed 21 July 2023. Retrieved from: <https://www.fda.gov/tobacco-products/market-and-distribute-tobacco-product/premarket-tobacco-product-applications>

Vidmantas U., Willeke, K., Grinshpun, S.A., Donnelly, J., Lin, X. & Mainelis, G. Aerosolization of Particles from a Bubbling Liquid: Characteristics and Generator Development, *AEROSOL SCI TECH*, 26:2, 175-190, DOI: 10.1080/02786829708965423, 1997.

405 Yahaya, A., Okoh, O. O., Agunbiade, F. O., & Okoh, A. I. Occurrence of phenolic derivatives in Buffalo River of Eastern Cape South Africa: Exposure risk evaluation. *ECOTOX ENVIRON SAFE*, 171, 887-893. doi:<https://doi.org/10.1016/j.ecoenv.2019.01.037>, 2019.

Zhao, T., Shu, S., Guo, Q., & Zhu, Y. Effects of design parameters and puff topography on heating coil temperature and mainstream aerosols in electronic cigarettes. *ATMOS ENVIRON*, 134. doi:10.1016/j.atmosenv.2016.03.027, 2016.

410

415

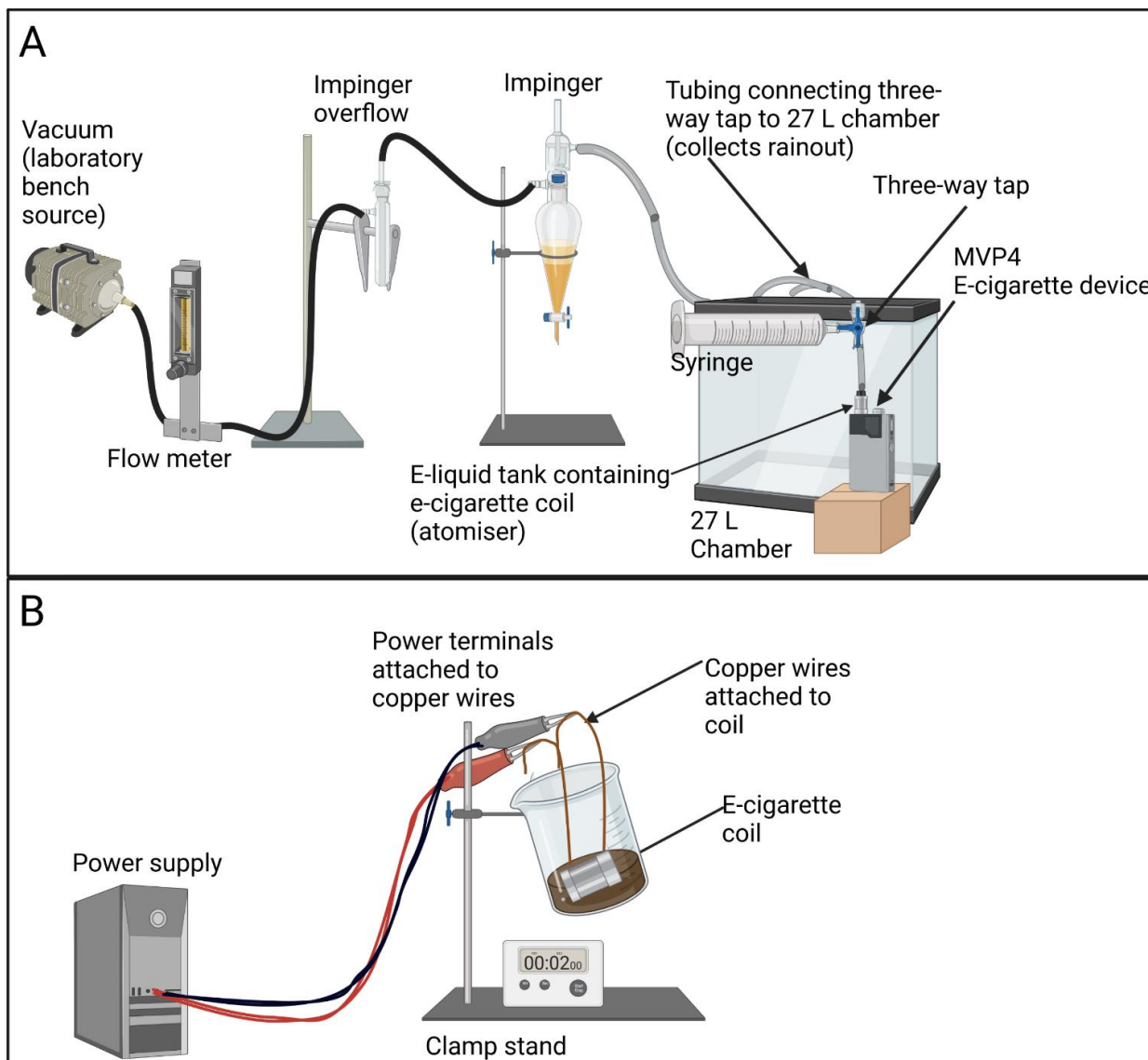


Figure 1. Vaping set-ups. A. Typical vaping set up. A vacuum drew air through the system at ~ 3 Lpm. The aerosol was drawn into a 60 mL syringe and a three-way tap was turned to allow the syringe to push the aerosol through two ~ 4 mm ID and ~ 15 cm tubing lengths and into the 27 L chamber for mixing. Air containing aerosol was drawn first into an impinger with 50:50 propylene glycol:glycerin base liquid. **B. Accelerated vaping set-up.** The power supply was attached to copper wires which were attached to an e-cigarette coil (Kanthal BVC) which was always submerged in e-liquid within a 100 mL beaker. Key differences between methods: 1) Aerosol is allowed to mix in air before capture rather than impinged immediately in liquid, 2) Volume of liquid is different: 3.5 mL is present in atomiser with typical method, compared to ~ 30 mL in beaker with accelerated method. Key similarities between methods: both methods apply the same amperage to the coil (and therefore heat coil to the same temperature).

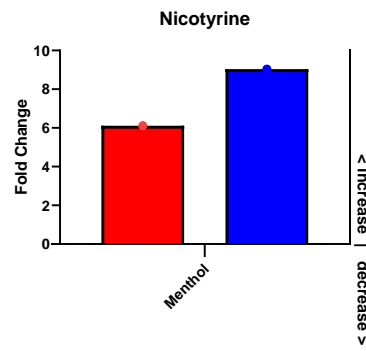
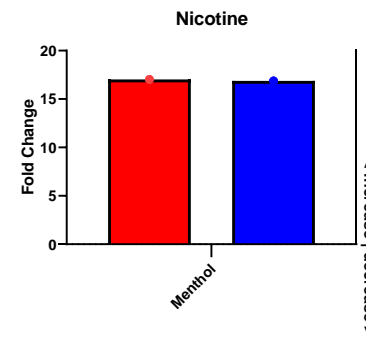
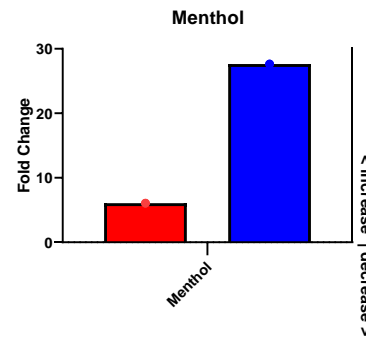
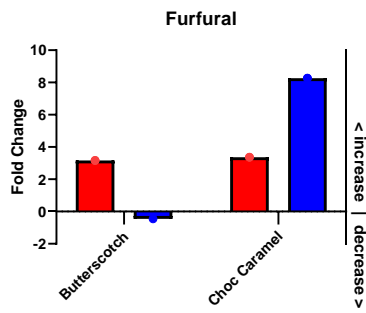
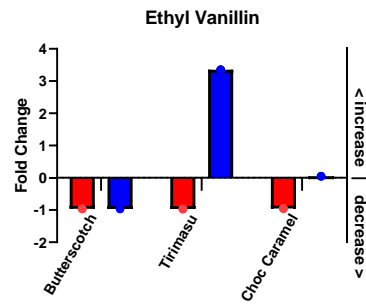
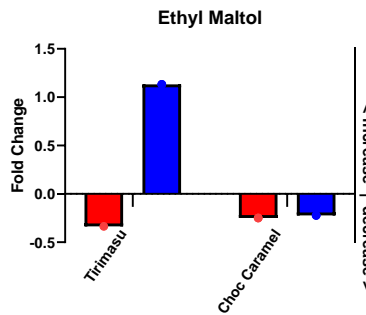
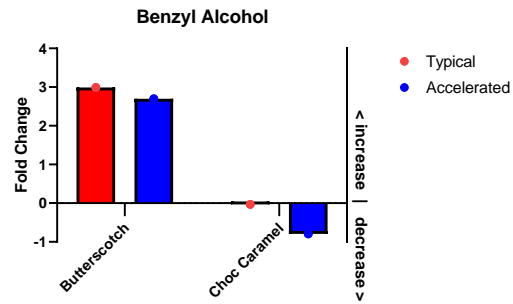
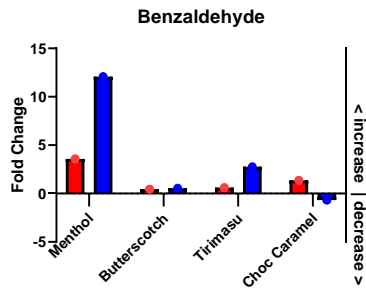


Figure 2. Fold change comparison between heating methods. E-liquid flavors are described on the X axis. Accelerated and Typical methods are indicated by the blue and red respectively. Y axis indicates the fold change compared to unheated i.e. Fold change = (Y-X)/X, where X is the unheated sample (mg/L concentration) and Y is the heated sample (mg/L concentration). Note different scales. Values that fall below the horizontal line at zero indicate a decrease in concentration from unheated sample, values above zero are an increase from unheated. If any chemical was not detected in a particular flavor then those flavors are not shown (e.g. there was no ethyl maltol detected in the menthol or butterscotch flavor).

Table 1. Assessment of 13 different chemical compounds from four e-liquids in both unheated and heated (accelerated and typical vaping methods). Unheated sample is comprised of only e-liquid in “fresh” or un-vaped form. Heated sample for typical method assessed the leftover from the atomiser tank and accelerated method assessed heated e-liquid produced from a coil submerged in e-liquid. Chemicals are listed in alphabetical order. Grey shading = increase from unheated sample. Acc = “accelerated” or “aged”. U = undetected. Italicized results indicate that the compound was not present in one of the heating method samples and therefore fold change analysis was not possible for that chemical in that e-liquid flavor. Numbers are accurate to three significant figures (digits) and values less than zero are displayed with logarithmic scale for ease of reading.

Chemical Tested (mg/L)	Menthol unheated	Menthol Heated (acc)	Menthol Heated (typical)	Butter-scotch unheated	Butter-scotch Heated (acc)	Butter-scotch Heated (typical)	Tiramisu unheated	Tiramisu heated (acc)	Tiramisu heated (typical)	Choc Caramel unheated	Choc Caramel heated (acc)	Choc Caramel heated (typical)
Benzaldehyde	9.87x10 ⁻³	1.29x10 ⁻¹	4.51x10 ⁻²	1.49	2.33	2.16	6.24x10 ⁻¹	2.48	1.06	1.49	2.81x10 ⁻¹	2.05
Benzyl Alcohol	u	u	u	13.40	49.53	53.46	345	u	285	238	49.53	231
Ethyl Maltol	u	u	u	u	u	u	339	723	225	891	693	669
Ethyl Vanillin	u	u	u	45.65	1.46	1.63	747	3250	23.73	1660	1740	76.98
Eugenol	u	u	u	u	u	u	<i>1.89x10⁻¹</i>	u	<i>1.25x10⁻¹</i>	u	u	u
Furfural	u	u	u	1.73x10 ⁻¹	9.38x10 ⁻²	7.21x10 ⁻¹	<i>6.57x10⁻¹</i>	u	<i>8.68x10⁻¹</i>	5.16x10 ⁻¹	4.78	2.25
Menthol	32.54	932	230	u	30.18	11.83	u	u	u	u	u	u
Nicotine	4.77x10 ⁻¹	8.53	8.60	u	u	u	u	u	u	u	u	u
Nicotyrine	2.70x10 ⁻⁴	2.71x10 ⁻³	1.92x10 ⁻³	u	u	u	u	u	u	u	u	u
Thymol	u	u	u	u	u	u	u	u	u	u	u	u
Trans-cinnamaldehyde	u	u	u	u	u	u	u	u	u	u	13.39	u
2-chloro-phenol	u	u	u	u	u	u	u	u	u	u	u	u
4-(4-methoxyphenyl)-2-butanone	u	u	u	u	u	u	u	u	u	u	u	u